

stryker

System 9 power tools

Built on trust. Empowering you.

Powered by the latest battery technology, our first-ever wireless charger and a handpiece designed to fully grasp your needs, System 9 is our most advanced power tool ever—for unmatched performance, whenever you're ready.



Just like its powerful legacy, System 9 power tools didn't start with us.

It started with you.

At Stryker, we've watched the mastery of a physician at work. We've stood inches away from millimeter precision. We've witnessed what the right tools make possible in the right hands—building trust time and time again. And yet a trust this deep can always be deeper, which is why we created System 9 power tools.

Every point of contact with System 9 power tools is intended to optimize ease of use and peace of mind. It includes our **first-ever wireless charger, advanced battery technology and an informed handpiece design that's been rethought top to bottom.**

Excellence isn't just about what happens when you pull the trigger. It's also about what happens before and beyond the procedure. From the **holistic ergonomic handpiece design powered through the innovative twist-lock connection to maintaining sterility and charge of batteries post-sterilization,** System 9 power tools enhance the entire experience—for the whole hospital.

Customer feedback is at the forefront of everything we do, while the voice of Dr. Homer Stryker is also in the back of our minds. That's how we unite the proven legacy of the past with the possibilities of the future—to build trust that empowers you.



**Sternum
saw**

**SABO
saw**

**Precision
oscillating
tip saw**

**Sagittal
saw**

**Reciprocating
saw**

System 9 power tools prove that we can enhance the most important power of all—the bond between a surgeon and the tool in their hand.



Cordless driver

Single trigger rotary drill

Dual trigger rotary drill

EZout powered acetabular revision system

Building trust that **empowers for 40 years**

Sometimes by tweaking, other times through total transformation, we've made our power tools better without changing the core of what made them the most dependable in the first place.

From the early days of Dr. Homer Stryker to the introduction of System 9, the same commitment has been clear. Together with our customers, we create products that address the current and future surgical landscape.



1983

OP 90
power
tools



1991

S2000
power
tools



1998

System 4
power
tools



2002

System 5
power
tools



**For four
decades,
we've
proactively
turned key
insights into
generation-
defining
enhancements**



2006

System 6

power
tools



2012

System 7

power
tools



2017

System 8

power
tools



2023

System 9

power
tools

Multiple handpieces, **all kinds of confidence**

A versatile selection to meet your procedural needs



Innovative technology **in your hands today**

Controlled, consistent and accurate

**EZout powered
acetabular revision
system**



**Precision oscillating
tip saw**



**Cordless driver with
POWEReam technology**



Cut with controlled power with the System 9 EZout system. The system powers the cutting, providing you with predictable, consistent results.

Enhanced design and expanded line of precision cartridges for ankle, shoulder, hip and knee arthroplasty, evolving the way you cut.

Drill, cut and ream with Stryker's exclusive technology, designed to provide a significant increase in torque for challenging procedures.

Now you have an **eye for power**

Efficiency through tactile and visual feedback

System 9's battery indicators are designed to be brighter and repositioned, creating greater visibility for surgeons and the surgical team. The battery housing features new contacts so that your battery is more dependable and durable.

Battery visibility, amplified

Indicator light upon connection




Solid green
Fully charged



Blinking amber
Low battery (<2 minute cut)



Solid amber
Fully depleted or error detected

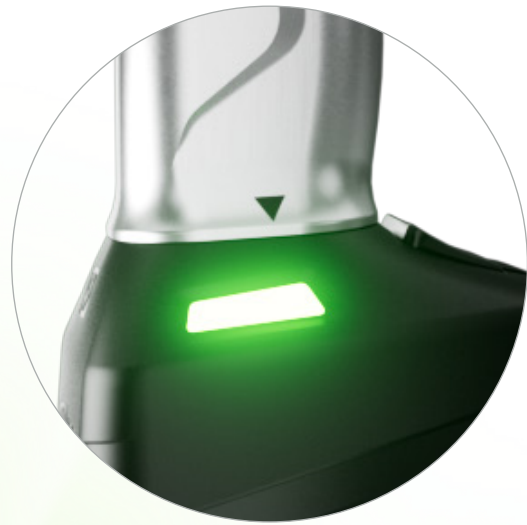

The charger that **changes everything**

Convenience meets competence

Our new wireless charger is the first of its kind—an advanced piece of technology for a smarter, more streamlined OR experience.

- ✔ Most versatile battery charging capabilities to date, allowing you to charge your batteries in any environment that accommodates your workflow
- ✔ Ability to confirm battery status pre-procedure
- ✔ Simplified user interface
- ✔ Ability to charge battery post sterilization while maintaining sterility, including ability to charge through blue wrap
- ✔ Custom top housing designed to accommodate various charging methods
- ✔ Foreign object detection application



A handle on design **and comfort**

**Built for confidence,
innovated for toughness**



Thoughtful contours

Well-balanced and designed to fit in your palm while reducing the number of pressure points while cutting.

.....

Built to outlast

Over time, exposure to pressurized water and harsh chemicals can wear down the toughest power tools—until now. Introducing our most robust anodized coating paired with IPX 9 rating for lasting performance.

.....

Improved safety pockets

Avoid accidental activation and perform with greater confidence and peace of mind.

.....

Twist-lock

Our new handpiece-to-battery interface provides tactile and visual feedback, enabling efficient assembly and confident connection.

We've got **you covered**

Cleaners and detergents

Proven to protect
your investment.

We've got a solution for all your needs, whether auto wash, manual wash, ultrasonic, point-of-use pretreatment or instrument lubrication.



Smart Equipment Management

Smart devices. Actionable data.
Informed decisions.

Stryker's Smart Equipment Management (SEM) is an online data module designed to help you optimize the management of your System 9 power tools through device-derived, actionable data.



ProCare Services

Protect your investment.

As the Original Equipment Manufacturer (OEM) we are experts in your new capital and your trusted partner offering:

- Proprietary product knowledge and original parts
- Expert technicians providing the ultimate customer support
- Loaner equipment to maximize uptime
- Flexible and unique service solutions for your business

Safety and reliability is at the heart of what we do, using an unauthorized service provider could put your patients at a greater risk and there is potential for a significant drop in performance and reliability of your new capital investment. Contact your Stryker Sales Rep today to get a customized ProCare Services plan for all of your Stryker equipment.



Flex Financial

Add flexibility to your budget.

Through Stryker's Flex Financial business, we provide a range of smart alternatives designed to fit your organization's needs. We offer flexibility beyond a cash purchase, enabling you to acquire our full portfolio of products. Ask your Sales Representative for assistance in finding the financial solution that best fits your needs.



Ordering information

Large bone power

9202-999-000	EZout System	3000-010-000	gallon (9.46L)	4100-120-000	T-Latch
9202-000-000	EZout Saw system		Auto62 Neutral Detergent;	4100-135-000	Hudson modified trinkle drill
9203-000-000	Single trigger rotary drill		low-foaming high-rinsing	4100-160-000	Trinkle drill
9205-000-000	Dual trigger rotary drill	3000-011-000	neutral concentrate	4100-210-000	Large AO reamer
9206-000-000	Reciprocating saw		2/2.5 gallon (9.46L)	4100-231-000	1/4" Reamer chuck with key
9207-000-000	Sternum saw		Auto62 Instrument	4100-235-000	Hudson modified
9208-000-000	Sagittal saw		Lubricant; neutral non-		trinkle reamer
9209-000-000	Stryker Precision Saw		silicone instrument	4100-260-000	Trinkle reamer
9215-000-000	System 9 battery pack, large		lubricant/2.5 gallon (9.46L)	4100-215-000	DHS/DCS reamer
9110-120-000US	System 9 charger			4100-355-000	Radiolucent right angle drive

Medium bone power

4605-000-000	Cordless driver	3000-003-000	Blu62 Pretreatment Foam -	4100-400-000	Sagittal saw attachment
4607-000-000	CD NXT handpiece	3000-007-000	10/32 oz (947mL)	4100-410-000	Right angle drive adaptor
4608-000-000	SABO Sagittal Saw	3000-002-000	Blu62 Pretreatment Foam	4100-600-000	1:1 Bur attachment
9212-000-000	System 9 battery pack, small		refill - 4/1 gallon (3.78L)	4100-700-000	High speed bur attachment
		3000-004-000	Blu62 Instrument detergent	4100-131-132	Universal chuck key
			- 20/8 oz (237mL)	1331-001-009	1/4" Chuck key
		3000-005-000	Blu62 Instrument detergent	1355-005-000	5/32" Chuck key
			is now 3/1 gallon (3.78L)		
		3000-001-000	Blu62 Instrument detergent		
			dispenser		
		3000-008-000	E62 Enzymatic detergent		
			- 4/1 gallon (3.78L)		

Aseptic

9126-100-000	System 9 non-sterile battery charger	3000-001-000	Blu62 Instrument detergent	8203-026-000	Dual trigger wire collet
9126-110-000	Large non-sterile battery		- 2/2.5 gallon (9.46L)		0.7 - 2.0mm
9222-110-000	Small non-sterile battery	3000-008-000	Blu62 Instrument detergent	8203-036-000	Single trigger wire collet
9126-120-000	Large aseptic battery housing		E62 Enzymatic detergent		0.7 - 2.0mm
9222-120-000	Small aseptic battery housing		- 4/1 gallon (3.78L)	8203-126-000	Dual trigger pin collet
9126-130-000	Large transfer shield				2.0 - 3.2mm
9222-130-000	Small transfer shield			8203-136-000	Single trigger pin collet
					2.0 - 3.2mm
				8203-226-000	Dual trigger large collet
					3.0 - 4.2mm
				8203-236-000	Single trigger large collet
					3.0 - 4.2mm
				8203-133-000	1/4" (6.4mm) Locking
					keyless chuck
				6203-110-000	Small AO quick connect
				6203-113-000	Hudson quick connect
				6203-131-000	1/4" Chuck with key
				6203-132-000	5/32" Chuck with key
				6203-135-000	Hudson modified trinkle quick connect
					Modified trinkle quick connect
				6203-150-000	connect
				6203-160-000	Trinkle quick connect
				6203-210-000	Large AO quick connect
				6203-215-000	DHS/DCS quick connect
				6203-170-000	Triathalon adaptor
				7207-003-000	Sternum blade guard
				7207-002-000	Sternum blade guard long

Smart Equipment Management

0280-100-010	SEM user license fee (qty 1)	4100-131-000	1/4" Drill chuck with key	6203-110-000	Small AO quick connect
0280-200-100	SEM handpieces - data subscription	4100-132-000	5/32" Drill chuck with key	6203-113-000	Hudson quick connect
0280-200-105	SEM batteries - data subscription	4405-210-000	POWEReam large AO	6203-131-000	1/4" Chuck with key
0280-300-120W	SEM System 9 charger charger - data subscription	4405-213-000	POWEReam Hudson	6203-132-000	5/32" Chuck with key
		4405-231-000	POWEReam 1/4" chuck with key	6203-135-000	Hudson modified trinkle quick connect
		4405-235-000	POWEReam Hudson modified trinkle		Modified trinkle quick connect
		4405-260-000	POWEReam trinkle		
		4100-126-000	Adjustable pin collet		
			2.0 - 3.2mm		
		4100-062-000	Wire collet 0.7 - 1.8mm		
		4100-125-000	Pin collet 2.0 - 3.2mm		
		4100-225-000	Pin collet 3.0 - 4.2mm		
		4100-110-000	Small AO drill		

Automatic cleaners and detergents

3000-009-000	Auto62 Enzymatic Cleaner; non-foaming high-rinsing enzymatic concentrate 2/2.5				
--------------	--	--	--	--	--

To order

Customer Service | P 800 253 3210 | F 877 790 7772 | inst.stryker.cs@stryker.com

Orthopaedic instruments

This document is intended solely for the use of healthcare professionals. A physician must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. We do not dispense medical advice and recommend that physicians be